Bruksanvisning / User Guide / Anleitung

Madrascover / Mattress Cover / Matratzenbezug ST

○ **MC 001-1915** [HMI: 35536]
ST - Small

Mått / Dimensions / Größe:
L 70 x B/W 205 cm

○ **MC 001-1925** [HMI: 35529]
ST - Medium

Mått / Dimensions / Größe:
L 135 x B/W 205 cm

○ **MC 001-1938** [HMI: 35530]
ST - Large

Mått / Dimensions / Größe:
L 200 x B/W 205 cm

○ **MC 001-2038** [HMI: 40022]
ST - X-Large

Mått / Dimensions / Größe:
L 200 x B/W 225 cm

PERNOVA

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Osäker?: +45 8693 8585

Washing- and tumble instructions

Washing temperature: **40 - 85°C**
Tumble dry: **max 60°C**
Airdrying: extends durability significantly

Do not use softener

Use detergent without bleach and chlorine-containing compounds

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UK

WARNING

There is a risk of the user slipping over the edge of the bed

Never leave the user lying on master **care's** products with lowered safety rails

Never leave the user with master **care's** products attached to the lift or with lowered safety rails

IMPORTANT

- These products must not be used if they show critical wear or are otherwise damaged
- Always use the products in accordance with this manual
- Always check the stitching on these products - both before and after use
- If in doubt about the proper use of the products, please contact master **care**/your dealer immediately
- master **care** assumes no liability for accidents resulting in injuries if these products are not used in accordance with this manual
- Trading with master **care** is conducted in accordance with the company's at all time applicable Trade and Business Conditions
- Velcro straps must be closed by washing and drying
- During wash and drying, the washing straps MUST be applied to the velcro on the top piece

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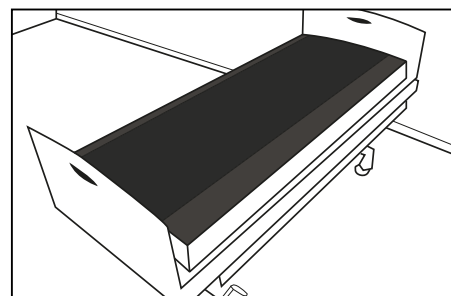
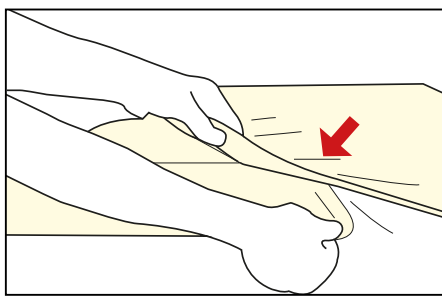
UK Place the mattress cover by the main body of the mattress. Hold the cover firmly to the mattress and fold the mattress upside down to the bed's foot

Then lock the ends of the mattress cover together with the attached Velcro fasteners and tighten. Now fold the mattress from the foot toward the main body and lock the ends in the same way



The mattress cover is now fitted and tightened. Checked continuously

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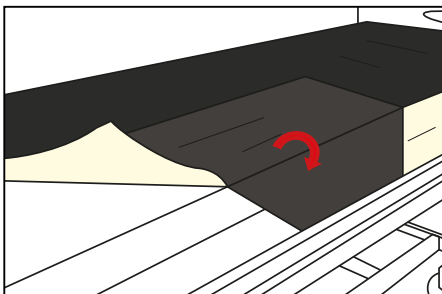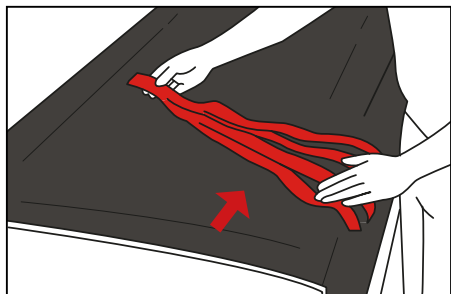UK If mattress covers are fitted on high mattresses, it is recommended that 4 mattress cover extensions are purchased

If a top mattress is placed on top of the ordinary mattress, fit the mattress cover around both mattresses

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UK

- Easy fitting
- 10 cm non-slip coating along the edge on both sides
- Unique sliding properties - HLF (Hyper Low Friction)
- Prevents pressure ulcers
- Extremely durable

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MC 001-1732 [HMI: 49368]
Master Turner Soft - **Large**

Mått / Dimensions / Größe:
L 215 x B/W 250

MC 001-2540 [HMI: 49370]
Master Turner Two In One Soft - **Large**

Mått / Dimensions / Größe:
L 215 x B/W 145 cm

MC 001-1740 [HMI: 84778]
Master Turner Pro Soft - **Large**

Mått / Dimensions / Größe:
L 215 x B/W 150 cm

MC 001-1700 [HMI: 49207]
Master Turner DUO Soft - Botten / Bottom /
Unterteil - **Large**

Mått / Dimensions / Größe:
L 160 x B/W 150 cm

MC 001-1701 [HMI: 49206]
Master Turner DUO Soft Topp / Top / Oberteil -
Large

Mått / Dimensions / Größe:
L 100 x B/W 150 cm


MC 001-1950 [HMI: 85234]
Master Turner Dragremmar / Locking Straps /
Schließschlaufen

Mått / Dimensions / Größe:
L 47 x B/W 4 cm

MC 001-2500 [HMI: 49873]
Master Turner DUO Two In One Soft Botten /
Bottom / Unterteil - **Large**

Mått / Dimensions / Größe:
L 160 x B/W 145 cm

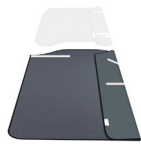MC 001-2501 [HMI: 49872]
Master Turner DUO Two In One Soft Topp / Top /
Oberteil - **Large**

Mått / Dimensions / Größe:
L 100 x B/W 145 cm

MC 001-1947 [HMI: 84373]
Master Turner Dragremmar / Pull Straps /
Zugband

Mått / Dimensions / Größe:
L 140 x B/W 5 cm

MC 001-8040 [HMI: 36382]
Master Comfort Stödkil / Support Wedge /
Stützkeile - 2 stk. / 2 pcs. / 2 Stück

Mått / Dimensions / Größe:
L 40 x B/W 15 x H 5 cm

MC 001-8050 [HMI: 34331]
Master Comfort Magkil - Kort /
Stomach Wedge - Short / Bauchkeil - Kurz

Mått / Dimensions / Größe:
L 43 x B/W 25 x H 14 cm

MC 001-8020 [HMI: 45939]
Master Comfort Tvättkil - Medium / Washing
Wedge / Waschkeil

Mått / Dimensions / Größe:
L 70 x B/W 40 x H 30 cm

MC 001-1948 [HMI: 81998]
Master Turner Justeringsremmar /
Adjustment Straps / Regulierungsurte

Mått / Dimensions / Größe:
L 40-80 x B/W 3 cm

MC 001-1203 [HMI: 34302]
Nylon Master Öppen - Medium / Open / Offen

Mått / Dimensions / Größe:
L 80 x B/W 160 cm

MC 001-1951 [HMI: 81134]
Mattress cover extender - **Large**

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Washing temperature: **40 - 85°C**

Tumble dry: **max 60°C**

Airdrying: extends durability significantly

Do not use softener.

Use detergent without bleach and chlorine-containing compounds.

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WARNING

There is a risk of the user slipping over the edge of the bed

Never leave the user lying on master **care's** products with lowered safety rails

Never leave the user with master **care's** products attached to the lift or with lowered safety rails

IMPORTANT

- These products must not be used if they show critical wear or are otherwise damaged
- Always use the products in accordance with this manual
- Always check the stitching on these products - both before and after use
- If in doubt about the proper use of the products, please contact master **care**/your dealer immediately

- master **care** assumes no liability for accidents resulting in injuries if these products are not used in accordance with this manual
- Trading with master **care** is conducted in accordance with the company's at all time applicable Trade and Business Conditions
- Velcro straps must be closed by washing and drying
- During wash and drying, the washing straps MUST be applied to the velcro on the top piece

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Pernova Hjälpmedel AB